

Stamina™

Tilt Chair and
Ophthalmic Stand

Unparalleled strength & style.



Metro 
Ophthalmic

Reichert®
TECHNOLOGIES

Advancing Eye Care. American Innovation. • www.reichert.com

AMETEK®

The Reichert® **Stamina™ Tilt Chair** meets the needs of the eye care professional who demands the very best in comfort, style, functionality, and durability.

Stamina™ Tilt Chair Features:




- **Premium quality.** Heavy-duty construction provides long-lasting performance.
- **Effortless maneuverability.** Counterbalanced 45° tilt and 360° rotation enables any staff member to easily position the patient as needed for different procedures and testing. Raise and lower the chair using the foot pedal, controls on the slit lamp arm or the stand control panel.
- **No more fluid leaks.** The industrial strength, Swiss designed linear actuator lift mechanism is quiet, durable and maintenance free, with a lift capacity of 560 lbs.
- **Easy access.** The stainless steel footrest and upholstered arms rotate upward and out of the way for convenient patient ingress and egress.
- **Maximum patient comfort.** Deep cushions offer a generous seat width, while premium upholstery remains soft to the touch without sacrificing durability.
- **Peace of mind.** 2 year standard warranty on the chair, 5 year warranty on linear actuator lift mechanism.

Stamina™ Tilt Chair Specifications

Catalog Number	15260	Rotation	360 degrees
Shipping Weight	295 lbs	Chair Tilt	45 degree
Gross Weight	260 lbs	Lower Limit	20.5 in.
Baseplate	22 in. W x 2.25 in. D	Upper Limit	30.5 in.
Lifting Column	560 lbs capacity	Input Voltage	120 VAC, 60 Hz
Patient Weight Limit	350 lbs capacity	Chair Fuse	8A (8 A 250V), 3AB, 3AG, RoHS
Lifting Column Stroke	11.8 in.		

The Reichert® **Stamina™ Ophthalmic Stand** meets the needs of the eye care professional who demands the very best in functionality, durability, and style.

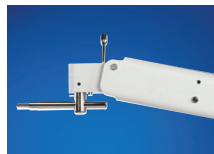

Stamina™ Ophthalmic Stand Features:

- **Strong, slim and stylish.** Durable construction outlasts the competition while contemporary styling complements your exam room. The slim profile and small footprint allows easy access to your patients and instruments without compromising floor space.
- **Logical control panel.** Buttons are located in close proximity to each other for easy access to the main on/off power switch, chair elevation and control, lower instrument arm, overhead lamp, projector and accessory instruments.
- **Convenience for accessories.** Three rechargeable instrument wells with charge indicators. Corded (BIO) instrument power and voltage selector. Three additional accessory outlets allow for consolidation of power cables.
- **Functionality at its finest.** The Phoroptor® arm comes fully counterbalanced for a Reichert® Ultramatic RX™ Phoroptor® and the lower instrument arm is counterbalanced for the Reichert® Xcel® 255 Slit Lamp. The built-in counterbalance adjustment mechanisms can be easily accessed in each arm. The height range on the lower arm accommodates wheelchairs and standard patient access.

Stamina™ Ophthalmic Stand Specifications

Catalog Number	15253
Height to top of pole:	68.5 in. (174 cm)
Base:	21.75 in. W x 22 in. D (52 cm x 56 cm)
Electrical	120 VAC, 50 - 60 Hz
Fuses	6 Amps Main
Corded Instrument Voltages	6 V, 7.5 V, 12 V, 12 V Fixed

Refractor Arm

Range in vertical motion:	12.5 in. (32 cm)
Length fully extended*:	26.25 in. (64 cm)

*Including Phoroptor® pole mount

Counterbalance weight limit: 20 lbs (9 kg)

Lower Arm

Range in vertical motion:	24 in. to 37.5 in. (61 cm x 95 cm)
Length fully extended:	35 in. (89 cm)
Counterbalance weight limit:	60 lbs (27 kg)

Optional Accessories

15254	Counterbalanced Keratometer Arm
15255	Non-counterbalanced Keratometer Arm
15258	Projector Mount
15002	Gas Lift Stool
15003	Gas Lift Stool with Backrest



Metro Ophthalmic

732-727-0607
www.metrophthalmic.com



Stamina™
Tilt Chair and Ophthalmic Stand

Unparalleled strength & style.



www.reichert.com

Call Reichert at +1 716-686-4500, toll-free 1-888-849-8955,
or contact your Authorized Reichert Technologies Distributor.
3362 Walden Avenue Depew, NY 14043 USA
© 2015 AMETEK, Inc. & Reichert, Inc. 15260-110-Rev A